

### 36-3502: Anti-Aurora B (Proliferation Marker) Monoclonal Antibody(Clone: AURKB/3121R)

<b>Clonality :</b>	Monoclonal
<b>Clone Name :</b>	AURKB/3121R
<b>Application :</b>	IHC
<b>Reactivity :</b>	Human
<b>Gene :</b>	AURKB
<b>Gene ID :</b>	9212
<b>Uniprot ID :</b>	Q96GD4
<b>Alternative Name :</b>	AIK2; AIM-1; ARK-2; AurB; AURKB; Aurora-1; Aurora and Ipl1 like midbody associated protein 1; Aurora kinase B; Aurora-B; Aurora-related kinase 2; Aurora/IPL1-related kinase 2; IPL1; Protein phosphatase 1 regulatory subunit 48 (PPP1R48); Serine/threonine-protein kinase 12; Serine/threonine-protein kinase aurora-B; STK1; STK12; STK5
<b>Isotype :</b>	Rabbit IgG
<b>Immunogen Information :</b>	Recombinant fragment (around aa 89-251) of human Aurora B protein (exact sequence is proprietary)

#### Description

Recognizes a protein of 39kDa, which is identified as Aurora B. The serine/threonine protein kinase aurora B (Aurora B) is a chromosomal passenger protein critical for accurate chromosome segregation, cytokinesis, protein localization to the centromere and kinetochore, correct microtubule-kinetochore attachment, and regulation of the mitotic checkpoint. Aurora B forms a tight complex with inner centrosome protein and survivin. Inactivation of any of these proteins causes similar defects in chromosome segregation. A significant overexpression of Aurora B has been found in a variety of human tumors including non-small cell lung carcinoma, astrocytoma, seminoma and carcinomas of the colon, prostate, endometrium and thyroid. The expression level of Aurora B is associated with cell proliferation and prognosis in these tumors.

#### Product Info

<b>Amount :</b>	20 µg / 100 µg
<b>Content :</b>	200 µg/ml of Ab Purified from Bioreactor Concentrate by Protein A/G. Prepared in 10mM PBS with 0.05% BSA & 0.05% azide. Also available WITHOUT BSA & azide at 1.0mg/ml.
<b>Storage condition :</b>	Antibody with azide - store at 2 to 8°C. Antibody without azide - store at -20 to -80°C. Antibody is stable for 24 months. Non-hazardous.

#### Application Note

Immunohistochemistry (Formalin-fixed) (1-2ug/ml for 30 minutes at RT)(Staining of formalin-fixed tissues requires heating tissue sections in 10mM Tris with 1mM EDTA, pH 9.0, for 45 min at 95&degC followed by cooling at RT for 20 minutes);

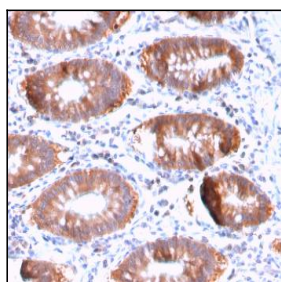


Fig. 1: Formalin-fixed, paraffin-embedded human Colon stained with Aurora B Rabbit Recombinant Monoclonal Antibody (AURKB/3121R).

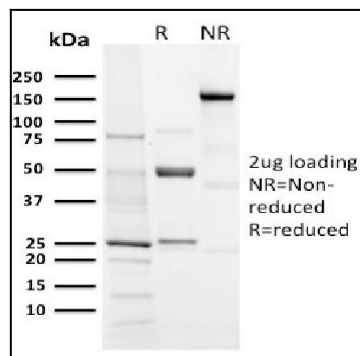


Fig. 2: SDS-PAGE Analysis Purified Aurora B Rabbit Recombinant Monoclonal (AURKB/3121R). Confirmation of Purity and Integrity of Antibody.