

9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982

Email: info@abeomics.com

## 36-1287: Monoclonal Antibody to Kappa Light Chain (B-Cell Marker)(Clone: HP6053)

Clonality: Monoclonal **Clone Name:** HP6053 FACS.WB.IHC Application: Reactivity: Human Gene: IGKV1D-16 **Uniprot ID:** P01601 Format: Purified **Alternative Name:** IGKV1D-16

**Isotype:** Mouse IgG1, kappa

Immunogen Information: Purified human Ig kappa chain

## **Description**

This MAb is specific to kappa light chain of immunoglobulin and shows no cross-reaction with lambda light chain or any of the five heavy chains. In mammals, the two light chains in an antibody are always identical, with only one type of light chain, kappa or lambda. The ratio of Kappa to Lambda is 70:30. However, with the occurrence of multiple myeloma or other B-cell malignancies this ratio is disturbed. Antibody to the kappa light chain is reportedly useful in the identification of leukemias, plasmacytomas, and certain non-Hodgkin's lymphomas. Demonstration of clonality in lymphoid infiltrates indicates that the infiltrate is malignant.

## **Product Info**

Amount:  $100 \mu g$ 

**Purification:** Affinity Chromatography

**Content:** 100 μg in 500 μl PBS containing 0.05% BSA and 0.05% sodium azide. Sodium azide is highly


toxic

**Storage condition :** Store the antibody at 4°C; stable for 6 months. For long-term storage; store at -20°C. Avoid

repeated freeze and thaw cycles.

## **Application Note**

Flow Cytometry (0.5-2.0ug/million cells); Western Blot (1-2ug/ml); Immunohistochemistry (Formalin-fixed) (1-2ug/ml for 30 min at RT)(Staining of formalin-fixed tissues requires heating tissue sections in 10mM Tris with 1mM EDTA, pH 9.0, for 45 min at 95&degC followed by cooling at RT for 20 minutes);

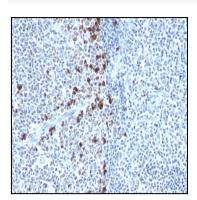

Formalin-fixed, paraffin-embedded human Tonsil stained with Kappa Light Chain Monoclonal Antibody (HP6053).