

9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982

Email: info@abeomics.com

## 36-1297: Monoclonal Antibody to Lambda Light Chain (B-Cell Marker)(Clone : Lamb14)(Discontinued)

Clonality: Monoclonal
Clone Name: Lamb14

**Application:** FACS,IF,WB,IHC

Reactivity: Human
Uniprot ID: P01701
Format: Purified
Isotype: Mouse IgG2a

Immunogen Information: Purified human lambda light chain

## **Description**

This MAb is specific to lambda light chain of immunoglobulin and shows no cross-reaction with lambda light chain or any of the five heavy chains. In mammals, the two light chains in an antibody are always identical, with only one type of light chain, kappa or lambda. The ratio of Kappa to Lambda is 70:30. However, with the occurrence of multiple myeloma or other B-cell malignancies this ratio is disturbed. Antibody to the lambda light chain is reportedly useful in the identification of leukemias, plasmacytomas, and certain non-Hodgkin's lymphomas. Demonstration of clonality in lymphoid infiltrates indicates that the infiltrate is malignant.

## **Product Info**

**Amount :** 100 μg

**Purification:** Affinity Chromatography

Content: 100 μg in 500 μl PBS containing 0.05% BSA and 0.05% sodium azide. Sodium azide is highly


toxic.

**Storage condition :** Store the antibody at 4°C; stable for 6 months. For long-term storage; store at -20°C. Avoid

repeated freeze and thaw cycles.

## **Application Note**

Flow Cytometry  $(0.5-1\tilde{A} \square \hat{A} \mu g/million$  cells in 0.1ml); Immunofluorescence  $(1-2\tilde{A} \square \hat{A} \mu g/ml)$ ; Western Blot  $(0.5-1\tilde{A} \square \hat{A} \mu g/ml)$ ; Immunohistology (Formalin-fixed)  $(0.5-1\tilde{A} \square \hat{A} \mu g/ml)$  for 30 minutes at RT); (Staining of formalin-fixed tissues requires boiling tissue sections in 10mM Citrate Buffer, pH 6.0, for 10-20 min followed by cooling at RT for 20 minutes); Optimal dilution for a specific application should be determined.

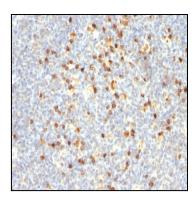


Formalin-fixed, paraffin-embedded human Tonsil stained with Lambda Light Chain Monoclonal Antibody (Lamb14).



9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982

Email: info@abeomics.com

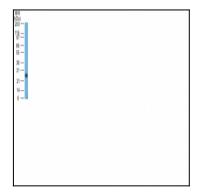

Western Blot of human Intestinal Lysate using Lambda Light Chain Monoclonal Antibody (Lamb14).