

9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982

Email: info@abeomics.com

36-1358: Monoclonal Antibody to Cytokeratin 8 (KRT8)(H1 + TS1)

Clonality: Monoclonal **Clone Name:** H1 + TS1 FACS.IF.IHC Application: Reactivity: Human Gene: KRT8 Gene ID: 3856 **Uniprot ID:** P05787 Format: Purified KRT8,CYK8 **Alternative Name:**

Isotype: Mouse IgG1, kappa + Mouse IgG1, kappa

Immunogen Information: Cytoskeleton preparation containing cytokeRatin 8 (H1); KeRatin preparation from a human

carcinoma (TS1)

Description

Cytokeratin 8 (CK8) belongs to the type II (or B or basic) subfamily of high molecular weight cytokeratins and exists in combination with cytokeratin 18 (CK18). CK8 is primarily found in the non-squamous epithelia and is present in majority of adenocarcinomas and ductal carcinomas. It is absent in squamous cell carcinomas. Hepatocellular carcinomas are defined by the use of antibodies that recognize only cytokeratin 8 and 18. CK8 exists on several types of normal and neoplastic epithelia, including many ductal and glandular epithelia such as colon, stomach, small intestine, trachea, and esophagus as well as in transitional epithelium. Anti-CK8 does not react with skeletal muscle or nerve cells. Epithelioid sarcoma, chordoma, and adamantinoma show strong positivity corresponding to that of simple epithelia (with antibodies against CK8, CK18 and CK19). Reportedly, anti-CK8 is useful for the differentiation of lobular (ring-like, perinuclear) from ductal (peripheral-predominant) carcinoma of the breast.

Product Info

Amount: $100 \mu g$

Purification: Affinity Chromatography

Content: 100 μg in 500 μl PBS containing 0.05% BSA and 0.05% sodium azide. Sodium azide is highly

toxic.

Storage condition : Store the antibody at 4°C; stable for 6 months. For long-term storage; store at -20°C. Avoid

repeated freeze and thaw cycles.


Application Note

Flow Cytometry (1-2ug/million cells); Immunofluorescence (1-2ug/ml); Immunohistochemistry (Formalin-fixed) (1-2ug/ml for 30 min at RT)(Staining of formalin-fixed tissues requires heating tissue sections in 10mM Tris with 1mM EDTA, pH 9.0, for 45 min at 95°C followed by cooling at RT for 20 minutes);



9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982

Email: info@abeomics.com

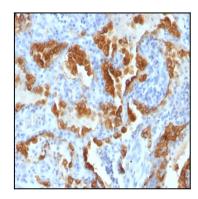

Formalin-fixed, paraffin-embedded Lung Carcinoma stained with Cytokeratin 8 Monoclonal Antibody (H1+TS1).