

## 12-1162: Anti-Cytokeratin 8/18 Recombinant Mouse Monoclonal Antibody (Clone:rKRT8.18/1346) (discontinued)

<b>Clonality :</b>	Monoclonal
<b>Clone Name :</b>	rKRT8.18/1346
<b>Application :</b>	IHC
<b>Reactivity :</b>	Human
<b>Gene :</b>	KRT8 & KRT18
<b>Gene ID :</b>	3856; 3875
<b>Uniprot ID :</b>	P05787; P05783
<b>Format :</b>	Purified
<b>Alternative Name :</b>	CARD2; CK8; CYK8; CYKER; Cytokeratin Endo A; DreK8; EndoA; K2C8; K8; Keratin 8; Krt 2.8; KRT8; Type-II Keratin Kb8. Cell Proliferation-inducing Gene 46 Protein; CK18; CYK18 Cytokeratin Endo B; K18; Keratin-18; Kerd; KRT18
<b>Isotype :</b>	Mouse IgG1, kappa
<b>Immunogen Information :</b>	Recombinant full-length human KRT8 protein and recombinant full-length human KRT18 protein

### Description

Cytokeratin 8 (CK8) belongs to the type II (or B or basic) subfamily of high molecular weight cytokeratins and exists in combination with cytokeratin 18 (CK18). This MAb cocktail recognizes all simple epithelia including glandular epithelium, for example thyroid, female breast, gastrointestinal tract, respiratory tract, and urogenital tract including transitional epithelium. All adenocarcinomas and most squamous carcinomas are positive but keratinizing squamous carcinomas are usually negative. This antibody is useful in demonstrating the presence of Paget cells; there is very little keratin 18 in the normal epidermis so only Paget cells are stained. Immunohistochemical staining with this MAb is indistinguishable from that obtained with monoclonal antibody 5D3.

### Product Info

<b>Amount :</b>	20 µg / 100 µg
<b>Purification :</b>	Protein A/G
<b>Content :</b>	200µg/ml of recombinant MAb purified by Protein A/G. Prepared in 10mM PBS with 0.05% BSA & 0.05% azide. Also available WITHOUT BSA & azide at 1.0mg/ml.
<b>Storage condition :</b>	Antibody with azide - store at 2 to 8°C. Antibody without azide - store at -20 to -80°C. Antibody is stable for 24 months. Non-hazardous.

### Application Note

Immunohistochemistry (Formalin-fixed) (1-2Âµg/ml for 30 minutes at RT),(Staining of formalin-fixed tissues requires heating tissue sections in 10mM Tris with 1mM EDTA, pH 9.0, for 45 min at 95Â°C followed by cooling at RT for 20 minutes)

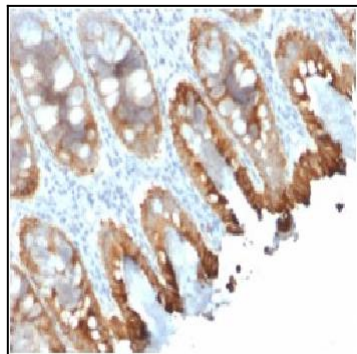


Figure 1: Formalin-fixed, paraffin-embedded Human colon stained with Cytokeratin 8/18 Mouse Recombinant Monoclonal Antibody (rKRT8.18/1346).

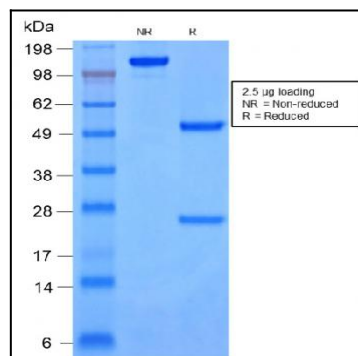


Figure 2: SDS-PAGE Analysis Purified Cytokeratin 8/18 Mouse Monoclonal Antibody (rKRT8.18/1346).