

14-513ACL: mCD73 Stable Cell Line-H

Application : Functional Assay

Description

mCD73 Stable Cell Line-H is a stably transfected HEK293 cell line which expresses mouse Cluster of Differentiation 73 (CD73 also known as ecto-5'-nucleotidase).

Sequence data: mCD73 (accession number NP_035981)

```
MRPAAAKVPKWLLLALSALLPQWPAASAWELTILHTNDVHSRLEQTSDDSTKCLNASLCV
GGVARLFTKVQQIRKEEPNVFLFDAGDQYQGTIWFTVYKGLEVAHFMMNILGYDAMALGNH
EFDNGVEGLIDPLLRNVKFPILSANIKARGPLAHQISGLFLPSKVLVSVGGEVVGIVGYTS
KETPFLSNPGTNLVFEDEISALQPEVDKLTNLVVKIILGHSGFEMDKLIAQKVRGVDI
VVGHSNTFLYTGNNPSKEVPAGKYPFIVTADDGRQVPVVQAYAFGKYLGYLKVEFDDKG
NVITSYGNPILLNSSIPEDATIKADINQWRIKLDNYSTQELGRTIVYLDGSTQTCRFREC
NMGNLICDAMINNNLRHPDEMFWNHVSMCIVNGGGIRSPIDEKNNGTITWENLAAVLPFG
GTFDLVQLKGSTLKKAFEHSVHRYGQSTGEFLQVGGIHVVYDINRKPWNRVVQLEVLCTK
CRVPIYEPLMDKVYKVTLPYLANGGDGFQMIKDELLKHDSGDQDISVVSEYISKMKVV
YPAVEGRIKFSAAASHYQGSFPLVILSWAMILILYQ
```

Product Info

Amount : For Profit / Non Profit
Content : Each vial contains 2 ~ 3 x 10⁶ cells in 1 ml of 90% FBS + 10% DMSO.
Storage condition : Immediately upon receipt, store in liquid nitrogen.

Application Note

Application:

- Screen for antibodies of mouse CD73 through Flow Cytometry.

Culture conditions:

Cells should be grown at 37°C with 5% CO₂ using DMEM medium (w/ L-Glutamine, 4.5g/L Glucose and Sodium Pyruvate) supplemented with 10% heat-inactivated FBS supplemented with 10% FBS and 1% Pen/Strep, plus 10 µg/ml of Blasticidin.

It is recommended to quickly thaw the frozen cells upon receipt or from liquid nitrogen in a 37°C water-bath, transfer to a tube containing 10 ml of growth medium without Blasticidin, spin down cells, resuspend cells in pre-warmed growth medium without Blasticidin, transfer resuspended cells to T25 flask and culture in 37°C-CO₂ incubator.

Leave the T25 flask in the incubator for 1~2 days without disturbing or changing the medium until cells completely recover viability and become adherent. Once cells are over 90% adherent, remove growth medium and passage the cells through trypsinization and centrifugation. At first passage, switch to growth medium containing Blasticidin. Cells should be split before they reach complete confluence.

To passage the cells, detach cells from culture vessel with Trypsin/EDTA, add complete growth medium and transfer to a tube, spin down cells, resuspend cells and seed appropriate aliquots of cells suspension into new culture vessels. Subcultivation ration = 1:10 to 1:20 weekly. To achieve satisfactory results, cells should not be passaged over 16 times.

LIMITED USE RESTRICTIONS:

THIS PRODUCT IS SOLELY FOR IN VITRO RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

By use of this product, user agrees to be bound by the terms of this limited use statement.

This product is solely for Internal Research Purposes and not for Commercial Purposes. Commercial Purposes include, but are not limited to (1) use of the cell line in manufacturing; (2) use of the cell line to provide a service, information or data; (3) use of the cell line for therapeutic, diagnostic or prophylactic purposes; or (4) resale of the cell line whether or not such cell lines are resold for use in research. The buyer cannot sell, give or otherwise transfer this product to a third party.

Commercial License Agreement is available for non-research use if applicable. Please contact Abeomics (info@abeomics.com).

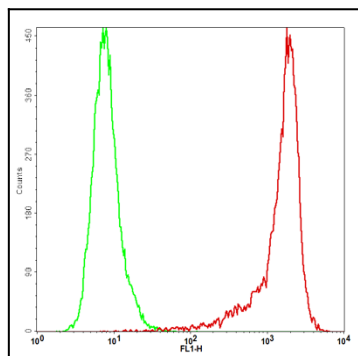


Fig-1: Detection of mouse CD73 in the HEK293/mCD73 stable cell line by Flow Cytometry [Cell surface staining]. HEK293 cells (Green); HEK293/mCD73 cells (Red).