

9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982

Email: info@abeomics.com

## 36-2772: Anti-Ki-67 (Proliferating Cell Marker) Monoclonal Antibody(Clone: MKI67/2462)

Clone Name: MKI67/2462

Application: FACS,IF,WB,IHC

 Reactivity :
 Human

 Gene :
 MKI67

 Gene ID :
 4288

 Uniprot ID :
 P46013

Alternative Name: KI-67; Ki67; KI-67 Antigen (KIA); MKI67; Proliferation related Ki-67 antigen

**Isotype:** Mouse IgG2b, kappa

Immunogen Information: Recombinant human Ki67 protein fragment (around aa 2293-2478) (exact seguence is proprietary)

## **Description**

Ki-67 antigen is a nuclear, non-histone protein that is present in all stages of the cell cycle except G0. This characteristic makes Ki-67 an excellent marker for proliferating cells and is commonly used as one of the prognostic factors in cancer studies. A correlation has been demonstrated between Ki-67 index and the histo-pathological grade of neoplasms. Assessment of Ki-67 expression in renal and ureter tumors shows a correlation between tumor proliferation and disease progression, thus making it possible to differentiate high-risk patients. Ki-67 expression may also prove to be important for distinguishing between malignant and benign peripheral nerve sheath tumors. Ki-67 labeling index has been shown to be a prognostic marker in a number of neoplasms including grade II astrocytoma, oligodendroglioma, colon carcinoma, and breast carcinoma. In general, Ki-67 is a good marker of proliferating cell populations.

## **Product Info**

**Amount :**  $20 \mu g / 100 \mu g$ 

Content: 200 µg/ml of Ab Purified from Bioreactor Concentrate by Protein A/G. Prepared in 10mM PBS with

0.05% BSA & 0.05% azide. Also available WITHOUT BSA & azide at 1.0mg/ml.

Storage condition:

Antibody with azide - store at 2 to 8°C. Antibody without azide - store at -20 to -80°C. Antibody is

stable for 24 months. Non-hazardous.

## **Application Note**

Flow Cytometry (1-2ug/million cells); Immunofluorescence (1-2ug/ml); Western Blot (1-2ug/ml); Immunohistochemistry (Formalin-fixed) (1-2ug/ml for 30 minutes at RT), (Staining of formalin-fixed tissues requires heating tissue sections in 10mM Tris with 1mM EDTA, pH 9.0, for 45 min at 95 &degC followed by cooling at RT for 20 minutes),

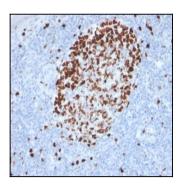


Fig. 1: Formalin-fixed, paraffin-embedded human Tonsil stained with Ki67-Monospecific Mouse Monoclonal Antibody (MKI67/2462).



Email: info@abeomics.com



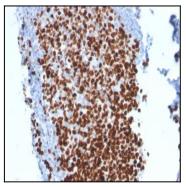


Fig. 2: Formalin-fixed, paraffin-embedded human Tonsil stained with Ki67-Monospecific Mouse Monoclonal Antibody (MKI67/2462).

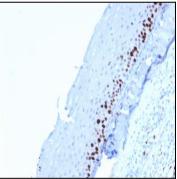


Fig. 3: Formalin-fixed, paraffin-embedded human Tonsil-Skin stained with Ki67-Monospecific Mouse Monoclonal Antibody (MKI67/2462).

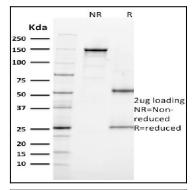


Fig. 4: SDS-PAGE Analysis Purified Ki67-Monospecific Mouse Monoclonal Antibody (MKI67/2462). Confirmation of Purity and Integrity of Antibody.

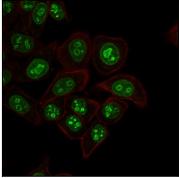


Fig. 5: Immunofluorescent staining of HeLa cells. Ki67-Monospecific Mouse Monoclonal Antibody (MKI67/2462) followed by goat anti-Mouse IgG-CF488 (green0. Cell membrane is labeled with Phalloidin (red).



9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982

Email: info@abeomics.com

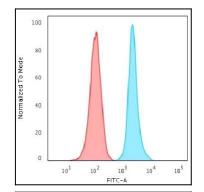


Fig. 6: Flow Cytometric Analysis of HeLa cells. Ki67-Monospecific Mouse Monoclonal Antibody (MKI67/2462) followed by goat anti-Mouse IgG-CF488 (blue); Isotype Control (red).

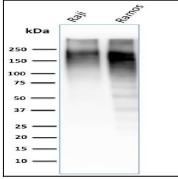


Fig. 7: Western Blot Analysis of Raji and Ramos cell lysate using Ki67-Monospecific Mouse Monoclonal Antibody (MKI67/2462).