∗ abeomics

36-2957: Anti-ACTH (Adrenocorticotrophic Hormone) (Pituitary Marker) Monoclonal Antibody(Clone: CLIP/1407)

Clonality :	Monoclonal
Clone Name :	CLIP/1407
Application :	IHC
Reactivity :	Human
Gene :	POMC
Gene ID :	5443
Uniprot ID :	P01189
Alternative Name :	Adrenocorticotropin; alpha or beta or gamma Melanocyte Stimulating Hormone (MSH) or Melanotropin; beta-Endorphin; beta or gamma Lipotropin (LPH); CLIP; Met Enkephalin; POC; POMC
Isotype :	Mouse IgG1, kappa
Immunogen Information : Synthetic peptide corresponding to aa25-39 of human ACTH (NGAEDESAEAFPLEF)	

Description

ACTH (same as Corticotropin) is a 39 amino acid active peptide produced by the anterior pituitary. This MAb is specific to CLIP (aa25-39 of ACTH); does not react with Synacthen (aa1-24 of ACTH). POMC (pro-opiomelanocortin or corticotropinlipotropin) is a 267 amino acid polypeptide hormone precursor that goes throµgh extensive, tissue-specific posttranslational processing by convertases. POMC is cleaved into ten hormone chains named NPP, ACTH, alpha-MSH (Melanocyte Stimulating Hormone), beta-MSH, gamma-MSH, CLIP (corticotropin-like intermediary peptide), Lipotropin-beta, Lipotropin-gamma, beta-endorphin and Met-enkephalin. ACTH is also produced by cells of immune system (T-cells, B-cells, and macrophages) in response to stimuli associated with stress. Anti-ACTH is a useful marker in classification of pituitary tumors and the study of pituitary disease. It reacts with ACTH-producing cells (corticotrophs).It also may react with other tumors (e.g. some small cell carcinomas of the lung) causing paraneoplastic syndromes by secreting ACTH.

Product Info

Amount :	20 μg / 100 μg
Content :	200 μ g/ml of Ab Purified from Bioreactor Concentrate by Protein A/G. Prepared in 10mM PBS with 0.05% BSA & 0.05% azide. Also available WITHOUT BSA & azide at 1.0mg/ml.
Storage condition :	Antibody with azide - store at 2 to 8°C. Antibody without azide - store at -20 to -80°C. Antibody is stable for 24 months. Non-hazardous.

Application Note

Immunohistochemistry (Formalin-fixed) (1-2ug/ml for 30 minutes at RT),(Staining of formalin-fixed tissues requires heating tissue sections in 10mM Tris with 1mM EDTA, pH 9.0, for 45 min at 95 & degC followed by cooling at RT for 20 minutes),

₩ abeomics

9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982 Email: info@abeomics.com

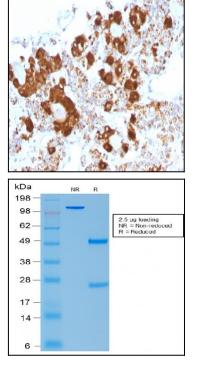


Fig. 1: Formalin-fixed, paraffin-embedded human Pituitary stained with ACTH Mouse Monoclonal Antibody (CLIP/1407).

Fig. 2: SDS-PAGE Analysis Purified ACTH Mouse Recombinant Monoclonal Antibody (rCLIP/1407). Confirmation of Purity and Integrity of Antibody.