

9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982

Email: info@abeomics.com

36-3624: Anti-CD44 / HCAM Std. Monoclonal Antibody(Clone: HCAM/2875R)

Clonality: Monoclonal
Clone Name: HCAM/2875R
Application: IHC,WB
Reactivity: Human
Gene: CD44
Gene ID: 960

LHR; BA-1; chondroitin sulfate proteoglycan 8 (CSPG8); Epican; Extracellular Matrix Receptor

III (ECM III); GP90 Lymphocyte Homing Adhesion Receptor; HCAM; HCELL; hematopoietic cell

Alternative Name : E- and L-selectin ligand; Heparan Sulfate Proteoglycan; Hermes Antigen; HAS; HUTCH I;

Hyaluronate Receptor; Indian blood group; Inlu Related p80 Glycoprotein; Ly 24; MDU2;

MDU3; MIC4; MUTCH I; Phagocytic Glycoprotein 1 (PGP-1)

Isotype: Rabbit IgG

Immunogen Information: Stimulated human leukocytes

P16070

Description

Uniprot ID:

Recognizes a cell surface glycoprotein of 80-95kDa (CD44) on lymphocytes, monocytes, and granulocytes (Leucocyte Typing Workshop V). Its epitope is resistant to digestion by trypsin and chymotrypsin. This MAb selectively interferes with lymphocyte binding to lymph node, mucosal and synovial endothelium. The CD44 family of glycoproteins exists in a number of variant isoforms, the most common being the standard 85-95kDa or hematopoietic variant (CD44s). Higher molecular weight isoforms are described in epithelial cells (CD44v), which are believed to function in intercellular adhesion and stromal binding. CD44 immunostaining is commonly used for the discrimination of urothelial transitional cell carcinoma in-situ from non-neoplastic changes in the urothelium.

Product Info

Amount: 20 μg / 100 μg

Content: 200 µg/ml of Ab Purified from Bioreactor Concentrate by Protein A/G. Prepared in 10mM PBS

with 0.05% BSA & 0.05% azide. Also available WITHOUT BSA & azide at 1.0mg/ml.

Storage condition : Antibody with azide - store at 2 to 8°C. Antibody without azide - store at -20 to -80°C. Antibody

is stable for 24 months. Non-hazardous.

Application Note

Western Blot (1-2ug/ml); Immunohistochemistry (Formalin-fixed) (0.25-0.5ug/ml for 30 minutes at RT)(Staining of formalin-fixed tissues requires heating tissue sections in 10mM Tris with 1mM EDTA, pH 9.0, for 45 min at 95°C followed by cooling at RT for 20 minutes);



9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982

Email: info@abeomics.com

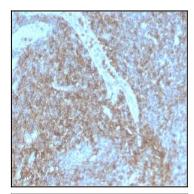


Fig. 1: Formalin-fixed, paraffin-embedded human Lymph Node stained with CD44 Rabbit Recombinant Monoclonal Antibody (HCAM/2875R).

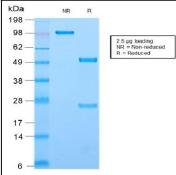


Fig. 2: SDS-PAGE Analysis Purified CD44 Rabbit Recombinant Monoclonal Antibody (HCAM/2875R). Confirmation of Integrity and Purity of Antibody.