

9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982

Email: info@abeomics.com

## 12-4119: Phospho-Stat5 (Tyr694) (Clone: B5) rabbit mAb

Clonality: Monoclonal Stat5Y694-B5 Clone Name:

Application: **FACS** 

Reactivity: Human, Mouse Conjugate: Unconjugated Format: Purified

Signal transducer and activator of transcription 5A, STAT5A, Signal transducer and activator of **Alternative Name:** 

transcription 5B, STAT5B

Rabbit IgG1k Isotype:

A synthetic phospho-peptide corresponding to residues surrounding Tyr694 of human phospho Immunogen Information:

## **Description**

Stat5 activation occurs in response to many ligands including prolactin, IL-2, growth hormone, and GM-CSF. Tyr694 phosphorylation is obligatory activation of Stat5 (1,2), and is mediated by Src upon erythropoietin stimulation (3). Phospho Stat5 is constitutively active in some leukemic cell types (4), and phospho Stat5 is found in some endothelial cells when treated with IL-3, suggesting its involvement in cell motility and angiogenesis (5). Stat5 has been shown to be encoded by two separate genes, Stat5a and Stat5b, which share over 90% amino acid sequence identity. In different cell types, Stat5a and Stat5b are independently regulated and activated. For example, interferon treatment predominantly activates Stat5a in U937 cells and Stat5b in HeLa cells (6).

## **Product Info**

Amount:  $20 \mu l / 200 \mu l$ 

Content: 1X PBS, 0.02% NaN3, 50% Glycerol, 0.1% BSA Store at -20°C. Avoid repeated freeze and thaw cycles. Storage condition:

## **Application Note**

1µg/mL - 0.001µg/mL. It is recommended that the reagent be titrated for optimal performance for each application. See product image legends for additional information.(0.5mg/ml)

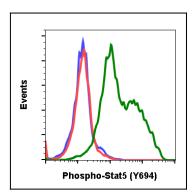


Fig-1: Flow cytometric analysis of K562 cells secondary antibody only negative control (blue) or treated with imatinib (red) or treated with pervanadate (green) using Phospho-Stat5 (Tyr694) antibody Stat5Y694-B5 at 0.05 µg/mL.



9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982

Email: info@abeomics.com

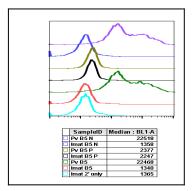


Fig 2 : Peptide blocking flow cytometric analysis of K562 cells secondary antibody only negative control (light blue) or treated with imatinib (red) or pervanadate-treated (green) or imatinib and blocked with phospho-peptide (black) or pervanadate and blocked with phospho peptide (gold) or imatinib and blocked with non-phospho peptide (dark blue) or pervanadate and blocked with non-phospho peptide (purple) using Phospho-Stat5 (Tyr694) antibody Stat5Y694-B5 at 0.1  $\mu$ g/mL.

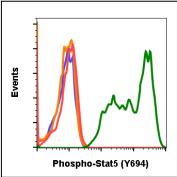


Fig-3: Flow cytometric analysis of NIH3T3 cells secondary antibody only negative control (blue) or 0.1  $\mu$ g/mL of isotype control (orange) or treated with imatinib (red) or with pervanadate (green) using Phospho-Stat5 (Tyr694) antibody Stat5Y694-B5 at 0.1  $\mu$ g/mL.

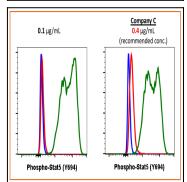


Fig-4: Flow cytometric analysis of 3T3 cells secondary antibody only negative control (blue) or treated with imatinib (red) or with pervanadate (green) using 0.1  $\mu$ g/mL of Phospho-Stat5 (Tyr694) antibody Stat5Y694-B5 or Company C antibody at 0.4  $\mu$ g/mL (manufacturer's recommended concentration).