

9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982

Email: info@abeomics.com

12-4344: Phospho-NFKB p65 (Ser468) (Clone: B9) rabbit mAb

Clonality: Monoclonal

Clone Name: NFKBp65S468-B9

Application :FACSReactivity :HumanConjugate :Unconjugated

Alternative Name: Transcription factor p65, Nuclear factor NF-kappa-B p65 subunit, Nuclear factor of kappa light

polypeptide gene enhancer in B-cells 3, RELA, NFKB3

Isotype: Rabbit IgG1k

Immunogen Information: A synthetic phospho-peptide corresponding to residues surrounding Ser536 of human

phospho-NFKB p65

Purified

Description

Format:

The nuclear factor kB (NFkB)/Rel family of transcription factors play a pivotal role in inflammatory and immune responses (1,2). NF-kappa-B is present in almost all cell types and is involved in many biological processes including immunity, inflammation, cell growth and differentiation, apoptosis, and tumorigenesis. NFkB is a homo- or heterodimeric complex formed by the Rel-like domain-containing proteins RELA/p65, RELB, NFkB1/p105, NFkB1/p50, REL and NFkB2/p52. The dimers bind at kB sites in the target gene DNA. Individual dimers have distinct preferences for different kB sites and can act as either transcriptional activators or repressors. NFkB Ser536 phosphorylation stimulates Lys310 acetylation and interaction of phospho NFkB with CBP. Acetylated/phospho NFkB induces enhanced transcriptional activity.

Product Info

Amount : $20 \mu l / 200 \mu l$

Content: 1X PBS, 0.09% NaN3, 0.2% BSA

Storage condition : Store at 2-8°C. Avoid repeated freeze and thaw cycles.

Application Note

For flow cytometric staining, the suggested use of this reagent is $5 \text{ A} \triangle \mu$ per million cells or $5 \text{ A} \triangle \mu$ per 100 A $\triangle \mu$ of staining volume. It is recommended that the reagent be titrated for optimal performance for each application. See product image legends for additional information. (0.5mg/ml, more than 200 western blots)

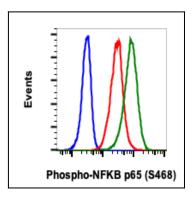


Fig-1: Flow cytometric analysis of A431 cells unstained and untreated with as negative control (blue) or untreated (red) or treated with staurosporine (green) and stained using phospho-NFKB p65 (Ser468) antibody NFKBP65S468-B9 at $0.01 \,\mu g/mL$.



9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982

Email: info@abeomics.com

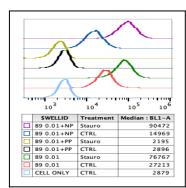


Fig 2 : Peptide blocking flow cytometric analysis of A431cells secondary antibody only negative control (light blue) or untreated (red) or treated with staurosporine (green) or untreated and blocked with phospho-peptide (black) or treated and blocked with phospho peptide (gold) or untreated and blocked with non-phospho peptide (dark blue) or treated and blocked with non-phospho peptide (purple) using Phospho-NFkB p65 (Ser468) antibody NFkBp65S468-B9 at 0.01µg/mL.