

9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982

Email: info@abeomics.com

30-2877AC: APC Conjugated Anti-Human CD303 Mab (Clone:15E3)

Clonality: Monoclonal

Clone Name: 15E3
Application: FACS
Reactivity: Human
Conjugate: APC
Gene: CD303
Gene ID: 170482
Uniprot ID: Q8WTT0

Alternative Name: BDCA2; CLEC4C; DLEC; HECL; CLECSF7; CLECSF11; PRO34150

Isotype: Mouse IgG1 kappa

Immunogen Information: CD303 ectodomain fused with human IgG Fc domain

Description

CD303 is an approximately 38 kDa type II transmembrane glycoprotein with an extracellular C-type lectin domain. It is a specific marker of plasmacytoid dendritic cells, and plays roles in capturing of pathogen-related oligosaccharide-containing antigens by them, and in their presentation to T cells. CD303 also mediates a potent inhibition of interferon alpha/beta production in plasmacytoid dendritic cells, thus it represents a potential target for lupus erythematosus therapy. Specificity:The mouse monoclonal antibody 15E3 recognizes an extracellular epitope of human CD303, a transmembrane glycoprotein expressed on plasmacytoid dendritic cells.

Product Info

Amount: 100 tests

Purified antibody is conjugated with activated allophycocyanin (APC) under optimum conditions

Purification: and unconjugated antibody and free fluorochrome are removed by size-exclusion

chromatography.

Content: Formulation: Stabilizing phosphate buffered saline (PBS), pH 7.4, 15 mM sodium azide

Storage condition: Store at 2-8°C. Protect from prolonged exposure to light. Do not freeze.

Application Note

Flow cytometry: The reagent is designed for analysis of human blood cells using 10 μ l reagent / 100 μ l of whole blood or 106 cells in a suspension. The content of a vial (1 ml) is sufficient for 100 tests.

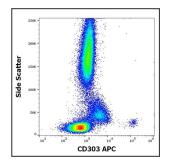


Figure 1: Flow cytometry surface staining pattern of human peripheral whole blood stained using anti-human CD303 (15E3) APC antibody (10 μ l reagent / 100 μ l of peripheral whole blood).



9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982

Email: info@abeomics.com

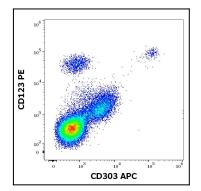


Figure 2: Flow cytometry multicolor surface staining pattern of human peripheral blood mononuclear cells stained using anti-human CD303 (15E3) APC antibody (10 μ l reagent / 100 μ l of peripheral whole blood) and anti-human CD123 (6H6) PE antibody (10 μ l reagent / 100 μ l of peripheral whole blood).

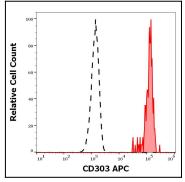


Figure 3: Separation of human CD303 positive CD123 positive plasmacytoid DC (red-filled) from neutrophil granulocytes (black-dashed) in flow cytometry analysis (surface staining) of human peripheral whole blood stained using anti-human CD303 (15E3) APC antibody (10 μ l reagent / 100 μ l of peripheral whole blood).