

30-2932PE: Streptavidin PE Conjugated

Application : FACS

Conjugate : PE

Description

Streptavidin PE is intended for detecting biotinylated antigens and biotinylated antibodies by FLOW cytometry.

Streptavidin PE is intended for detecting biotinylated antigens and biotinylated antibodies by flow cytometry.

Product Info

Amount :	100 ug
Purification :	Purified recombinant protein is conjugated with R-phycoerythrin (PE) under optimum conditions. Unconjugated protein and free fluorochrome are removed by size-exclusion chromatography.
Content :	Streptavidin PE - 1 vial (0.5 ml) of Ready-to-use solution. No reconstitution is necessary. Reagent is provided in stabilizing phosphate buffered saline (PBS) solution containing 15mM sodium azide.
Storage condition :	Store at 2-8°C. Protect from prolonged exposure to light. Do not freeze.

Application Note

Flow cytometry: Recommended use: 1 µg/ml. Extracellular and intracellular staining.

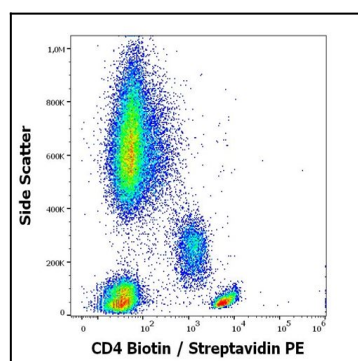


Figure 1: Flow cytometry surface staining pattern of human peripheral whole blood stained using anti-human CD4 (MEM-241) Biotin antibody and secondary stained using Streptavidin PE (concentration in sample 1 µg/ml).

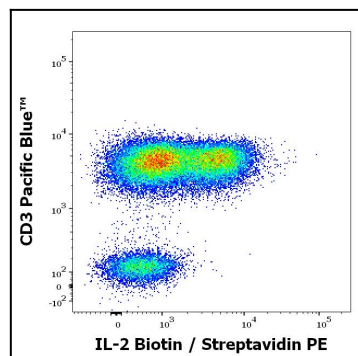


Figure 2: Flow cytometry multicolor surface staining pattern of PMA + Ionomycin stimulated and Brefeldin A treated human lymphocytes using anti-human CD3 (UCHT1) Pacific Blue™ antibody (4 µl reagent / 100 µl of peripheral whole blood) and intracellular staining of human lymphocytes using anti-human IL-2 (35C3) Biotin antibody (concentration in sample 0.5 µg/ml) and secondary stained using Streptavidin PE (concentration in sample 1 µg/ml).

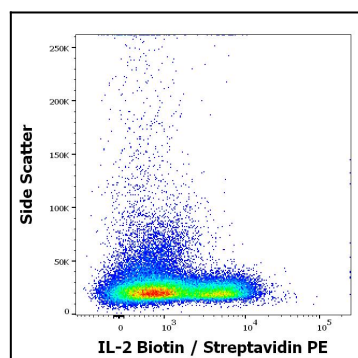


Figure 3: Flow cytometry intracellular staining pattern of PMA + Ionomycin stimulated and Brefeldin A treated human peripheral whole blood stained using anti-human IL-2 (35C3) Biotin antibody and secondary stained using Streptavidin PE (concentration in sample 1 µg/ml).