

9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982

Email: info@abeomics.com

## 36-1980: Monoclonal Antibody to CD79a (B-Cell Marker)(Clone: HM47/A9)

Clonality: Monoclonal Clone Name: HM47/A9

**Application :** FACS,IF,WB,IHC **Reactivity :** Human, Mouse, Rat

 Gene :
 CD79A

 Gene ID :
 973

 Uniprot ID :
 P11912

 Format :
 Purified

Alternative Name : CD79A,IGA,MB1
Isotype : Mouse IgG1, kappa

Immunogen Information: A synthetic peptide corresponding to aa 202-216 (GTYQDVGSLNIADVQ) of human CD79a

protein.

## **Description**

A disulphide-linked heterodimer, consisting of mb-1 (or CD79a) and B29 (or CD79b) polypeptides, is non-covalently associated with membrane-bound immunoglobulins on B cells. This complex of mb-1 and B29 polypeptides and immunoglobulin constitute the B cell Ag receptor. CD79a first appears at pre B cell stage, early in maturation, and persists until the plasma cell stage where it is found as an intracellular component. CD79a is found in the majority of acute leukemias of precursor B cell type, in B cell lines, B cell lymphomas, and in some myelomas. It is not present in myeloid or T cell lines. Anti-CD79a is generally used to complement anti-CD20 especially for mature B-cell lymphomas after treatment with Rituximab (anti-CD20). This antibody will stain many of the same lymphomas as anti-CD20, but also is more likely to stain B-lymphoblastic lymphoma/leukemia than is anti-CD20. Anti-CD79a also stains more cases of plasma cell myeloma and occasionally some types of endothelial cells as well.

## **Product Info**

**Amount :** 100 μg

**Purification:** Affinity Chromatography

Content: 100 μg in 500 μl PBS containing 0.05% BSA and 0.05% sodium azide. Sodium azide is highly

toxic.

Storage condition:

Storage condition:

Store the antibody at 4°C; stable for 6 months. For long-term storage; store at -20°C. Avoid

repeated freeze and thaw cycles.


## **Application Note**

Flow Cytometry (1-2ug/million cells); Immunofluorescence (1-2ug/ml); Western Blot (1-2ug/ml);Immunohistochemistry (Formalin-fixed) (1-2ug/ml for 30 minutes at RT)(Staining of formalin-fixed tissues requires heating tissue sections in 10mM Tris with 1mM EDTA, pH 9.0, for 45 min at 95&degC followed by cooling at RT for 20 minutes);



9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982

Email: info@abeomics.com

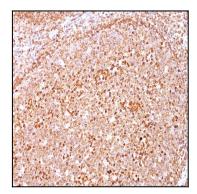

Formalin-fixed, paraffin-embedded human Tonsil stained with CD79a Ab (HM47/A9).