

9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982

Email: info@abeomics.com

## 36-11042: Monoclonal Antibody to Golgi Complex (Marker for Human Cells)(Clone: SPM581)

Clonality: Monoclonal Clone Name: SPM581

**Application:** FACS,IF,WB,ICC,IHC

Reactivity: Human Format: Purified

**Isotype:** Mouse IgG1, kappa

Immunogen Information: SU-DHL-1 large cell lymphoma cells.

## **Description**

This MAb recognizes Golgi complex in human cells. It is a part of a new panel of reagents, which recognizes subcellular organelles or compartments of human cells. These markers may be useful in identification of these organelles in cells, tissues, and biochemical preparations. It recognizes an antigen associated with the Golgi complex in human cells only. It can be used to stain the Golgi complex in cell or tissue preparations and can be used as a Golgi marker in subcellular fractions. It produces a diffuse staining pattern of the Golgi zone in normal and malignant cells and may be used to stain Golgi complex of cells in frozen tissue sections. It can also be used with paraformaldehyde fixed frozen tissue or cell preparations. This MAb is an excellent marker for human cells in xenographic model research. It reacts specifically with human cells.

## **Product Info**

Amount:  $100 \mu g$ 

**Purification:** Affinity Chromatography

**Content:** 100 μg in 500 μl PBS containing 0.05% BSA and 0.05% sodium azide. Sodium azide is highly


toxic.

**Storage condition :** Store the antibody at 4°C; stable for 6 months. For long-term storage; store at -20°C. Avoid

repeated freeze and thaw cycles.

## **Application Note**

Flow Cytometry (0.5-1.0ug/million cells);Immunofluorescence (1-2ug/ml); Western Blot (1-2ug/ml);Immunocytochemistry (Acetone or paraformaldehyde fixed) (1-2ug/ml for 30 min);Immunohistochemistry (Formalin-fixed) (1-2ug/ml for 30 minutes at RT)(Staining of formalin-fixed tissues requires heating tissue sections in 10mM Tris with 1mM EDTA, pH 9.0, for 45 min at 95&degC followed by cooling at RT for 20 minutes);

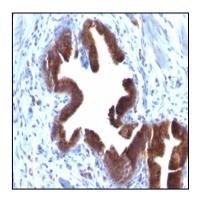

Formalin-fixed, paraffin-embedded human Gallbladder stained with Golgi Monoclonal Antibody (SPM581).