

9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982

Email: info@abeomics.com

36-11084: Monoclonal Antibody to Cytokeratin, pan (Epithelial Marker)(SPM115 + SPM116)

Clonality: Monoclonal

Clone Name: SPM115 + SPM116 Application: FACS.IF.WB.IHC Reactivity: Human, Mouse, Rat

Gene: KRT77 Gene ID: 374454 **Uniprot ID:** Q7Z794 Format: Purified KRT77,KRT1B **Alternative Name:** Mouse IgG1, kappa Isotype:

Immunogen Information: Human epidermal keRatin

Description

Twenty human keratins are resolved with two-dimensional gel electrophoresis into acidic (pl 6.0) subfamilies. This antibody cocktail recognizes acidic (Type I or LMW) and basic (Type II or HMW) cytokeratins, which 67kDa (CK1); 64kDa (CK3); 59kDa (CK4); 58kDa (CK5); 56kDa (CK6); 52kDa (CK8); 56.5kDa (CK10); 50kDa (CK14); 50kDa (CK15); 48kDa (CK16); 40kDa (CK19). Many studies have shown the usefulness of keratins as markers in cancer research and tumor diagnosis. AE-1/AE-3 is a broad spectrum anti pan-cytokeratin antibody cocktail, which differentiates epithelial tumors from non-epithelial tumors e.g. squamous vs. adenocarcinoma of the lung, liver carcinoma, breast cancer, and esophageal cancer. It has been used to characterize the source of various neoplasms and to study the distribution of cytokeratin containing cells in epithelia during normal development and during the development of epithelial neoplasms. This antibody stains cytokeratins present in normal and abnormal human tissues and has shown high sensitivity in the recognition of epithelial cells and carcinomas.

Product Info

Amount: 100 µg

Purification: Affinity Chromatography

100 µg in 500 µl PBS containing 0.05% BSA and 0.05% sodium azide. Sodium azide is highly Content:

toxic.

Store the antibody at 4°C; stable for 6 months. For long-term storage; store at -20°C. Avoid Storage condition :

repeated freeze and thaw cycles.


Application Note

Flow Cytometry (1-2ug/million cells); Immunofluorescence (1-2ug/ml); Western Blot (1-2ug/ml for 2 hours at RT); Immunohistochemistry (Formalin-fixed) (1-2ug/ml for 30 min at RT)(Staining of formalin-fixed tissues requires heating tissue sections in 10mM Tris with 1mM EDTA, pH 9.0, for 45 min at 95°C followed by cooling at RT for 20 minutes);



9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982

Email: info@abeomics.com

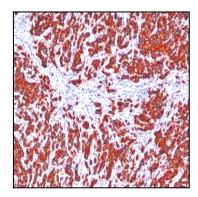

Formalin-fixed, paraffin-embedded human Breast Carcinoma stained with pan Cytokeratin Monoclonal Antibody cocktail (SPM115 \pm SPM116).