

## 14-132ACL: SRE LEEPORTER™ Luciferase Reporter-HEK293 Cell Line

**Application :** Functional Assay

### Description

The SRE LEEPORTER™ Luciferase Reporter cell line is a stably transfected HEK 293 cell line which expresses Renilla luciferase reporter gene under the transcriptional control of the serum response element (SRE). The SRE reporter cell line is designed to monitor MAPK/ERK activity and can be used for studying GPCR-linked MAPK/ERK signaling pathways as well as screening of agonists, antagonists or signaling inhibitors related with the MAPK/ERK signaling pathways. Functional activity of the cell line has been validated by serum stimulation assay (Figure 1).

### Product Info

**Amount :** 1 Vial  
**Content :** Each vial contains  $2 \sim 3 \times 10^6$  cells in 1 ml of 90% FBS + 10% DMSO.  
**Storage condition :** Immediately upon receipt, store in liquid nitrogen.

### Application Note

#### Potential application:

- Monitor the GPCR-linked MAPK/ERK signaling pathway.
- Screen for activators or inhibitors of the GPCR-linked MAPK/ERK signaling pathway.
- (Note: Abeomics did not test the potential applications above in-house.)

#### Culture conditions:

Cells should be grown at 37°C with 5% CO<sub>2</sub> using DMEM medium (w/ L-Glutamine, 4.5g/L Glucose and Sodium Pyruvate) supplemented with 10% heat-inactivated FBS and 1% Pen/Strep, plus 3 µg/ml of Puromycin (Note: Puromycin can be omitted during the reporter cell assays).

It is recommended to quickly thaw the frozen cells upon receipt or from liquid nitrogen in a 37°C water-bath, transfer to a tube containing 10 ml of growth medium without Puromycin, spin down cells, resuspend cells in pre-warmed growth medium without Puromycin, transfer resuspended cells to T25 flask and culture in 37°C-CO<sub>2</sub> incubator.

Leave the T25 flask in the incubator for 1~3 days without disturbing or changing the medium until cells completely recover viability and become adherent. Once cells are over 90% adherent, remove growth medium and passage the cells through trypsinization and centrifugation. At first passage, switch to growth medium containing Puromycin. Cells should be split before they reach complete confluence.

To passage the cells, detach cells from culture vessel with Trypsin/EDTA, add complete growth medium and transfer to a tube, spin down cells, resuspend cells and seed appropriate aliquots of cells suspension into new culture vessels. Subcultivation ration = 1:10 to 1:20 weekly. To achieve satisfactory results, cells should not be passaged over 16 times.

#### Functional validation:

A. Response of SRE LEEPORTER™ - HEK293 cells to fetal bovine serum (FBS).

1. Harvest SRE Leeporter™ – HEK293 cells and seed cells into a white solid-bottom 96-well microplate at  $5 \times 10^4$  cells/well in 100  $\mu$ l of growth medium containing 0.5% FBS to start serum starvation.
2. Incubate cells at 37°C in a CO<sub>2</sub> incubator for overnight.
3. The next day, stimulate cells with 40% FBS or 40% FBS plus phorbol 12-myristate 13-acetate (PMA).
4. Incubate at 37°C in a CO<sub>2</sub> incubator for 16 hours.
5. Equilibrate the plate to room temperature for 10 minutes.
6. Add 50  $\mu$ l of luciferase assay reagent (Abeomics, Cat #17-1101; Refer to the reagent datasheet for the detailed luciferase assay protocol) per well.
7. Read the plate in 1-5 minutes to measure luminescence using a microplate luminometer.

#### LIMITED USE RESTRICTIONS:

**THIS PRODUCT IS SOLELY FOR IN VITRO RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.**

**By use of this product, user agrees to be bound by the terms of this limited use statement.**

**This product is solely for Internal Research Purposes and not for Commercial Purposes. Commercial Purposes include, but are not limited to (1) use of the cell line in manufacturing; (2) use of the cell line to provide a service, information or data; (3) use of the cell line for therapeutic, diagnostic or prophylactic purposes; or (4) resale of the cell line whether or not such cell lines are resold for use in research. The buyer cannot sell, give or otherwise transfer this product to a third party.**

**Commercial License Agreement is available for non-research use if applicable. Please contact Abeomics ([info@abeomics.com](mailto:info@abeomics.com)).**

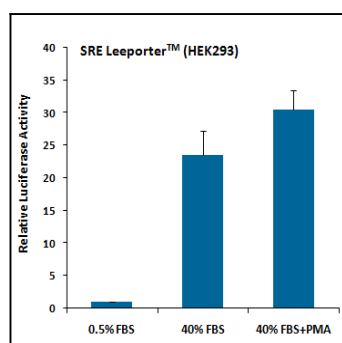


Fig-1: Induction of SRE activity by FBS or FBS plus PMA in SRE Leeporter™ - HEK293 cells.